

BioSpectrum GENIUS MORCELLATOR

Fast and precise



The Simai Genius Morcellator supplied by BioSpectrum Ltd has proven to be a “game changer” in terms of speed and efficiency. With key features including:

- Lightweight ergonomic handle design
- Adjustable blade speed to return between 13-20 grams per minute
- Fast consistent suction ensuring tissue retention throughout
- Single foot pedal offering 2 stage operation – suction, suction and morcellation
- High quality and durable cutting blades available in single use or reusable option
- Large capacity collection canisters up to 8 liters
- 12L Single Use Closed Fluid System also available

BioSpectrum Ltd

Unit 9 Halifax Court,
Fernwood Business Park,
Newark, Nottinghamshire,
NG24 3JP

Contact us for more information

<http://www.bio-spectrum.co.uk/>
 01636 705 084
customers@bio-spectrum.co.uk
https://twitter.com/Bio_Spectrum
<https://www.linkedin.com/company/bio-spectrumuk>
<https://www.youtube.com/channel/UCBtwuYGwGDZld8EA-hWLSXA>

Follow us



BioSpectrum GENIUS MORCELLATOR

Fast and precise



KEY SPECIFICATION

Blade length	360mm length, 4.8mm diameter
Blade options	Single Use or Reusable
Suction pressure	0-700 mmGh (0-90 Kpa)
Collection Cannisters	2 x 4 litre collection cannisters
Tissue collection	Yes contained in cannister 1

With a performance of between 13-20 grams per minute **the Genius Morcellator** really stands out from the competition. Well-built with a dedicated trolley the system is designed to be compact and elegant offering trouble free performance.

Arrange an evaluation

01636 705084
customers@bio-spectrum.co.uk



Fast: 13-20 g/m



Single Use or reusable blades



Precise



Constant suction



Ergonomic

BioSpectrum Ltd

Unit 9 Halifax Court,
Fernwood Business Park,
Newark, Nottinghamshire,
NG24 3JP

Contact us for more information

<http://www.bio-spectrum.co.uk/>
 01636 705 084
customers@bio-spectrum.co.uk
https://twitter.com/Bio_Spectrum
<https://www.linkedin.com/company/bio-spectrumuk>
<https://www.youtube.com/channel/UCBtwuYGwGDZld8EA-hWLSXA>

Follow us




BioSpectrum HoLEP Instrument Set

01. HD CYSTOSCOPE




- Ø 2.9 mm, autoclavable, 302 mm WL
- 30° / 12° viewing angle
- Probe channel: 1.2 mm
- Code: AA-003-038 / AA-003-044

02. WORKING ELEMENT

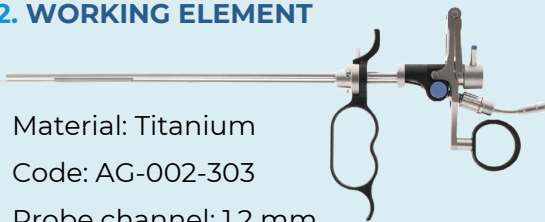

- Material: Titanium
- Code: AG-002-303
- Probe channel: 1.2 mm
- Passive cutting, handle closed, 19/22 Fr.

03. OUTER SHEATH

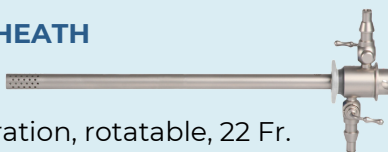

- With perforation, rotatable, 22 Fr.
- Code: AG-002-302

04. INNER SHEATH



- With ceramic isolation, rotatable
- Code: AG-002-301

05. OBTURATOR




- Standard, 19 Fr.
- Code: AA-003-011

06. VIEWING OBTURATOR




- 19 Fr.
- Code: AA-003-012

07. NEPHROSCOPE

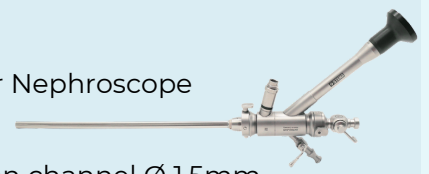

- 45° angled ocular Nephroscope
- 250mm length
- Ø 7.45mm, suction channel Ø 1.5mm, working channel Ø 4.7mm
- Code: AA-006-031

08. OUTER SHEATH (SECOND TYPE)

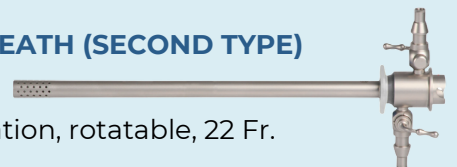

- With perforation, rotatable, 22 Fr.
- Code: AG-002-302

09. ACCESSORIES

- Adapter for laser probes → AG-100-006
- Touny-Borst Adapter → AG-002-073
- Sealing Cap for Nephroscope → AA-003-044

BioSpectrum Ltd

Unit 9 Halifax Court,
Fernwood Business Park,
Newark, Nottinghamshire,
NG24 3JP

Contact us for more information

<http://www.bio-spectrum.co.uk/>
 01636 705 084
customers@bio-spectrum.co.uk
https://twitter.com/Bio_Spectrum
<https://www.linkedin.com/company/bio-spectrumuk>
<https://www.youtube.com/channel/UCBtwuYGwGDZld8EA-hWLSXA>

Follow us





BioSpectrum

----- Creating a difference in healthcare -----

CONTACT US



01636 705 084



www.bio-spectrum.co.uk



customers@bio-spectrum.co.uk