

YOU ASKED FOR IT
YOU WAITED FOR IT

NOW, YOU HAVE IT!



PL-1000

HIGH RESOLUTION IMAGE

- Video processor panel equipped with image adjustment buttons to adapt or capture images

Light weight, easy to carry around

- No additional monitor required
- Real-time image provided via picture-in-picture during freeze



PL-2100

PL-1000				PL-2100			
Power	Voltage	100-240 V		Voltage	100-240 V		
	Frequency	50/60 Hz		Frequency	50/60 Hz		
	Rated input	45 VA		Rated input	48 VA		
Illumination	Lamp	2W LED		Lamp	2WLED		
	Lamp life	30,000 hrs		Lamp life	30,000 hrs		
	Color temperature	3.000K-7.000K		Color temperature	3.000K-7.000K		
Dimensions (W*D*H)	350mm*425mm*130mm,8kg			270mm*72mm*260mm,3kg			
Function	White balance,Sharpness,Freeze			White balance,Sharpness,Freeze			
Output	DVK1,S-Video*1,CVBS*1			HDMI*2			
Storage	USB 2.0			USB 2.0			
Screen	/			10.1 inch			


Accessories

Disposable Biopsy Forceps



Model	Working diameter	Length
FTL-BF-10-250-O-C-R	1.0mm,oval	2500mm
FTL-BF-10-290-O-C-R	1.0mm,oval	2900mm
FTL-BF-15-250-O-C-R	1.5mm,oval	2500mm
FTL-BF-15-290-O-C-R	1.5mm,oval	2900mm

Disposable Extraction Basket

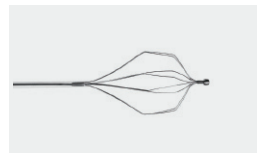

Model	Working diameter	Length
FTL-GF-10-250-15-FAA	1.0mm,opening range 15mm	2500mm
FTL-GF-10-290-15-FAA	1.0mm,opening range 15mm	2900mm

Disposable Snare




Model	Working diameter	Length
FTL-GF-10-290-9-O	1.0mm,opening range 25mm	2900mm

Disposable Distal End Tape Hood



Model	L*W*H	Color
WDA-160-20-01	160*20*0.1	Transparent
WDA-160-20-02	160*20*0.2	Transparent
WDA-160-20-03	160*20*0.3	Transparent
WDA-160-20-02B	160*20*0.2	Black
WDA-160-20-03B	160*20*0.3	Black

Disposable Endoscope Suction Button Set

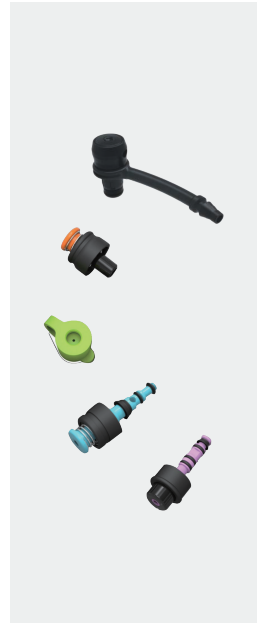

Model	Description	Specification
SG-MAJ-207	Suction valve(button)	Suitable for Olympus BF, ENF, CHF, CYF type slim endoscope
SG-DVK-OL-200/ SG-DVK-OL-200(W)	Olympus A/W Suction Button Set	For Olympus endoscopes (240,260,290series) 1.Air/water button; 2.Suction button; 3.Biopsy valve; 4.(Adapter)
SG-DVK-FJ-700/ SG-DVK-FJ-700(W) SG-DVK-FJ-500 SG-DVK-FJ-500(W)	Fuji A/W Suction Button Set	For Fuji endoscopes (601,700series) 1.Air/water button; 2.Suction button; 3.Biopsy valve; 4.(Adapter)
SG-DVK-PT-11090/ SG-DVK-PT-11090(W)	Pentax A/W Suction Button Set	For Pentax endoscopes (i10.etc. series) 1.Air/water button; 2.Suction button; 3.Biopsy valve; 4.(Adapter)



BRiViEW™

DISPOSABLE CHOLANGIOSCOPE
- for ERCP use


Shanghai SeeGen Photoelectric Technology Co., Ltd.

SEEGEN
vision borderless



Shanghai SeeGen Photoelectric Technology Co.,Ltd
Add.:Building 5,No.4277 Yindu Rd.Minhang Dist., 201108, Shanghai China
www.seegenmed.com
info@seegen.com.cn +86 021 5429 5355

SEEGEN
vision borderless

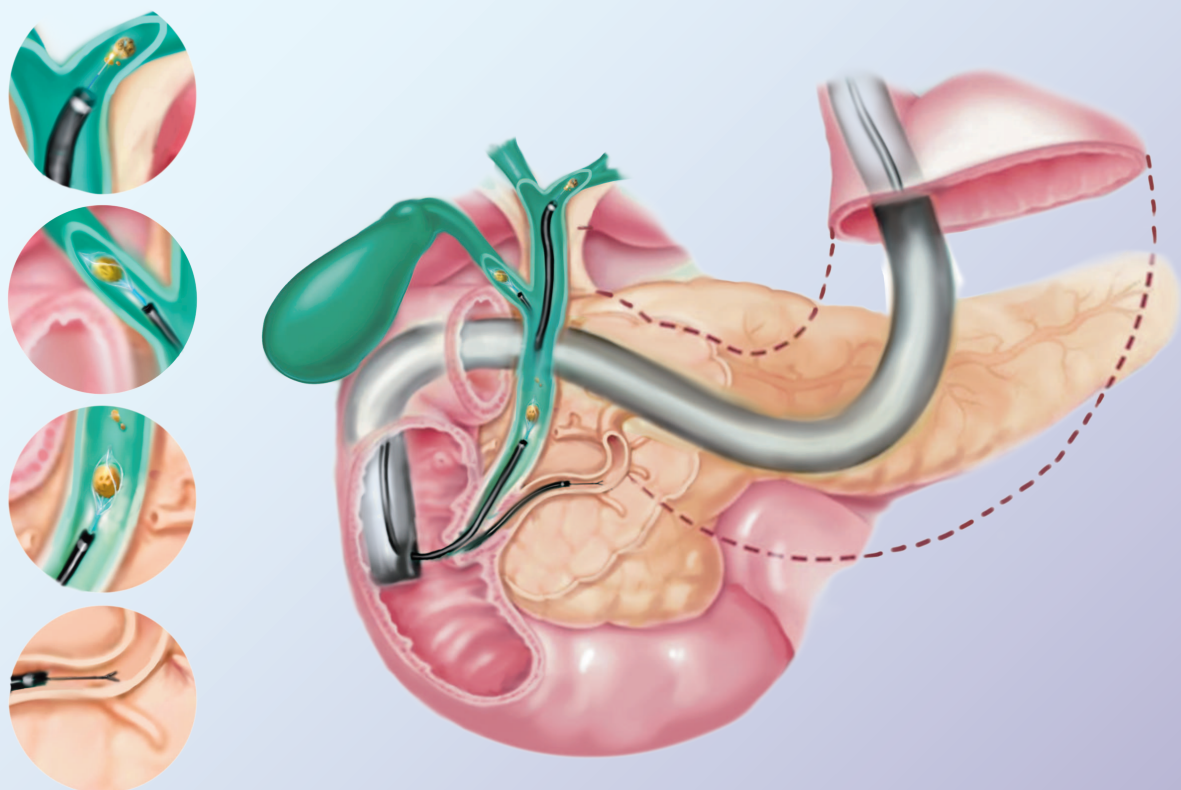

Total Solution For Biliary Diagnosis and Treatment



BRIVIEW DISPOSABLE CHOLANGIOSCOPE for ERCP use

Remote buttons for
image information obtainment


Instrument

4-way deflection for treatment
and complex cases

Working with Duodenoscope



Models	CHV-US100J-U	CHV-US110J-U	CHV-US120J-U	CHV-110J-U
Working length	2000mm	2000mm	2150mm	2000mm
Insertion tube width	3.3mm	3.5mm	2.6mm	3.8mm
Working channel diameter	1.2mm	1.2mm	0.9mm	1.8mm
Angle of deflection	up/down≥120°	up/down/left/right≥120°	up/down/left/right≥120°	up/down/left/right≥120°
Field of view	120°	120°	120°	120°
Depth of field	5-50mm	5-50mm	5-50mm	5-50mm
Working channel direction				

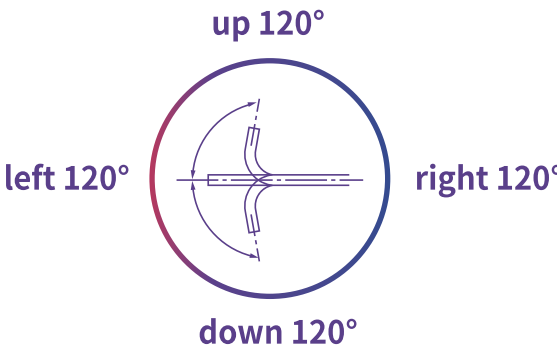


Water jet

Irrigation

Flow control for enhanced visibility

Suction

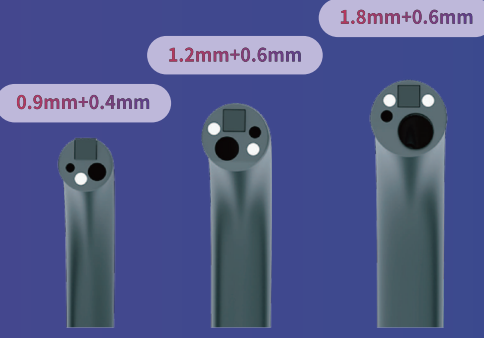

MULTI-SIZED WORKING CHANNELS

- 0.9mm slim channel: Slim design compatible with smaller duodenoscopes
Easy pancreatic duct, CBD and hepatic duct access
Tissue damages minimization
- 1.2mm working channel: For daily use
- 1.8mm working channel: Maximized specimen decreases biopsy frequency
Comprehensive accessory selection for cost efficiency